


AL #	Description	Size
AL-62364	Gel-Foam Cushion	16" x 16" x 2.5"
AL-62365	Gel-Foam Cushion	18" x 16" x 2.5"
AL-62366	Gel-Foam Cushion	20" x 16" x 2.5"
AL-62367	Thin-Line Gel Foam Cushion	18" x 16" x 1.5"

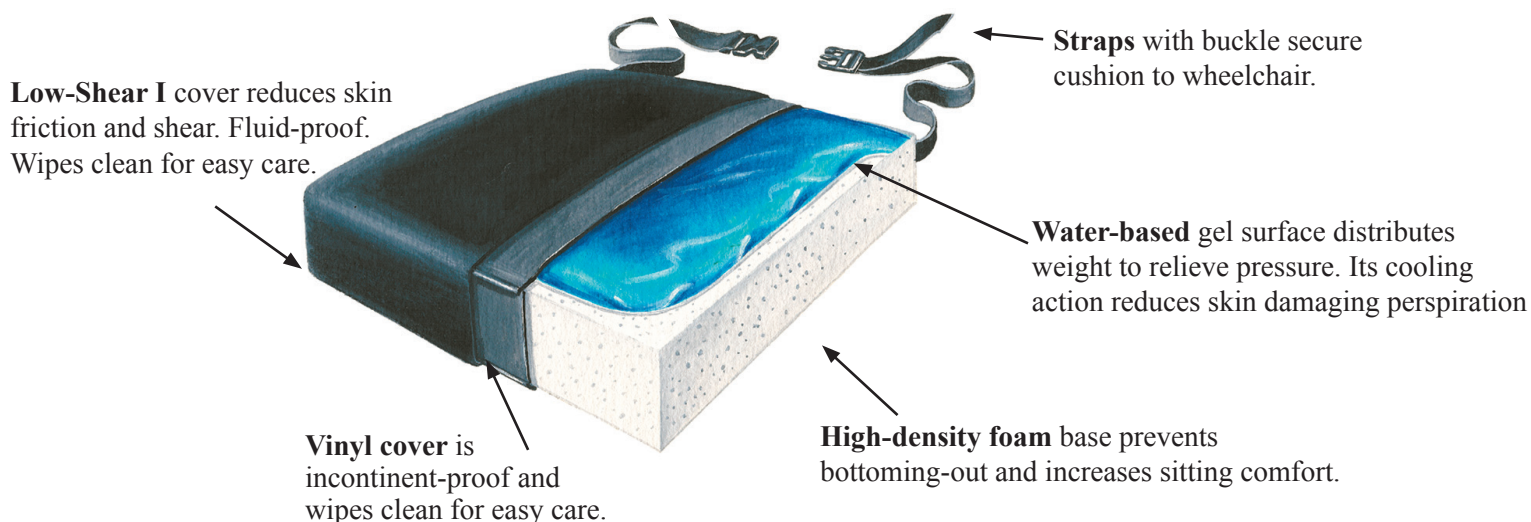

### Purpose

Gel-Foam Cushion is a combination gel and foam cushion designed to relieve pressure. The water based gel surface provides cooling action that reduces skin-damaging perspiration.

### Important Safety Reminder

The Low-Shear cover has a safety strap that prevents the cushion from slipping off the seat. This strap must be secured to the wheelchair frame.

### Easy Care

The vinyl cover may be wiped clean using a damp cloth and spray-type cleanser. A mild disinfectant may be used. The Low-Shear outer cover can be laundered and dried at temperatures that do not exceed 180°F (82°C). Do not use bleach.