


Conform Cushion  
AL# 62377



AL #	Description	Size
AL-62377	Conform Cushion w/Low Shear II Cover	18" x 16" x 1.5"

**Low-Shear** fluid proof cover reduces skin friction and shear and wipes clean for easy care.

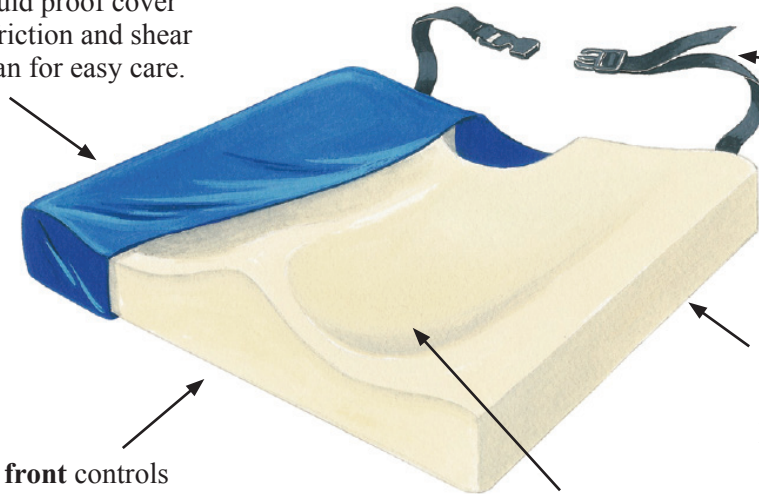

**Straps** with buckle secure cushion to wheelchair.

**Raised front** controls forward sliding.

**Temperature-sensitive visco-foam** conforms to user for maximum weight distribution and pressure reduction.

**Leg troughs** promote thigh alignment for postural symmetry



Visco-Elastic Foam

### Purpose

Made with temperature-sensitive visco-foam, this cushion provides pressure relief, supports appropriate posture and controls forward sliding.

### Important Safety Reminder

The Conform Cushion with a cloth outer cover has straps with buckles that must be secured to the wheelchair to prevent the cushion from sliding off the seat.

### Easy Care

The fluid proof cover may be wiped clean using a damp cloth and spray-type cleanser. A mild disinfectant may be used. Do not use bleach.

*This product is guaranteed for 6 months.*

Phone: 800-323-4282  
Fax: 800-950-1167



Like us on Facebook at:  
facebook.com/ALCOSalesServiceCo



Follow us on Twitter at:  
twitter.com/ALCOSales

info@alcosales.com  
www.alcosales.com