


Low-Shear I cover reduces skin friction and shear. Fluid-proof. Wipes clean for easy care.

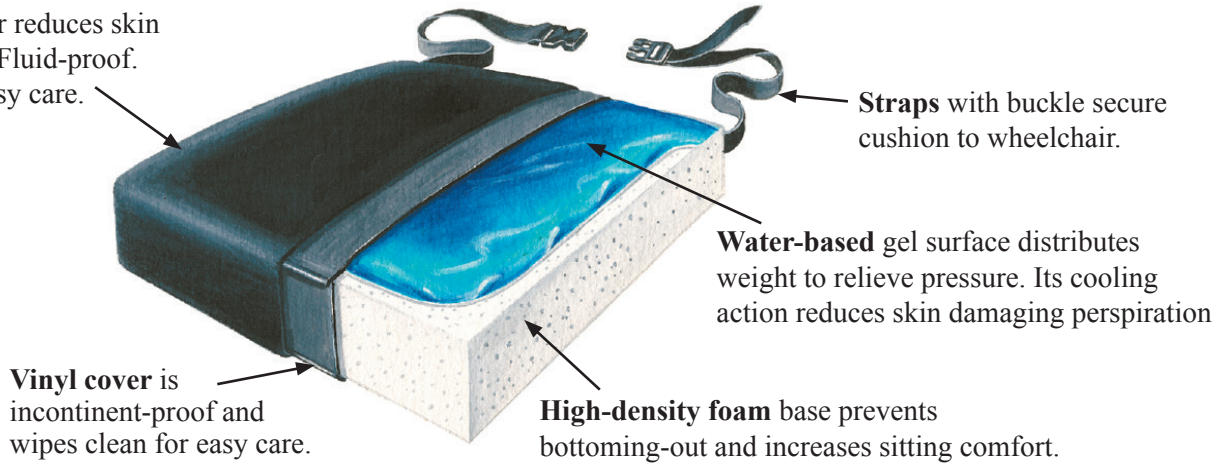

Straps with buckle secure cushion to wheelchair.

Water-based gel surface distributes weight to relieve pressure. Its cooling action reduces skin damaging perspiration

Vinyl cover is incontinent-proof and wipes clean for easy care.

High-density foam base prevents bottoming-out and increases sitting comfort.

AL #	Description	Size	Unit
AL-84698-16	Gel-Foam Thin-Line Vinyl Cushion w/ LSI Cover	16"W x 16"D x 1.5"H	EA
AL-84698-18	Gel-Foam Thin-Line Vinyl Cushion w/ LSI Cover	18"W x 16"D x 1.5"H	EA
AL-84698-20	Gel-Foam Thin-Line Vinyl Cushion w/ LSI Cover	20"W x 16"D x 1.5"H	EA

Purpose

The Gel-Foam Cushion is a gel over foam cushion designed to relieve pressure and increase comfort. The water based gel surface provides cooling action that reduces skin-damaging perspiration.

Easy to Use

Place cushion on a level surface and secure to the chair with safety ties.

Easy Care

The vinyl cover may be wiped clean using a damp cloth and spray-type cleanser. A mild disinfectant may be used. The Low-Shear outer cover can be laundered and dried at temperatures up to 180°F (82°C). Do not use bleach.